


Compact Scale

# HJ-150

*HJ-150 compact scale is suitable to carry with you anywhere in your pocket, bag or briefcase!*

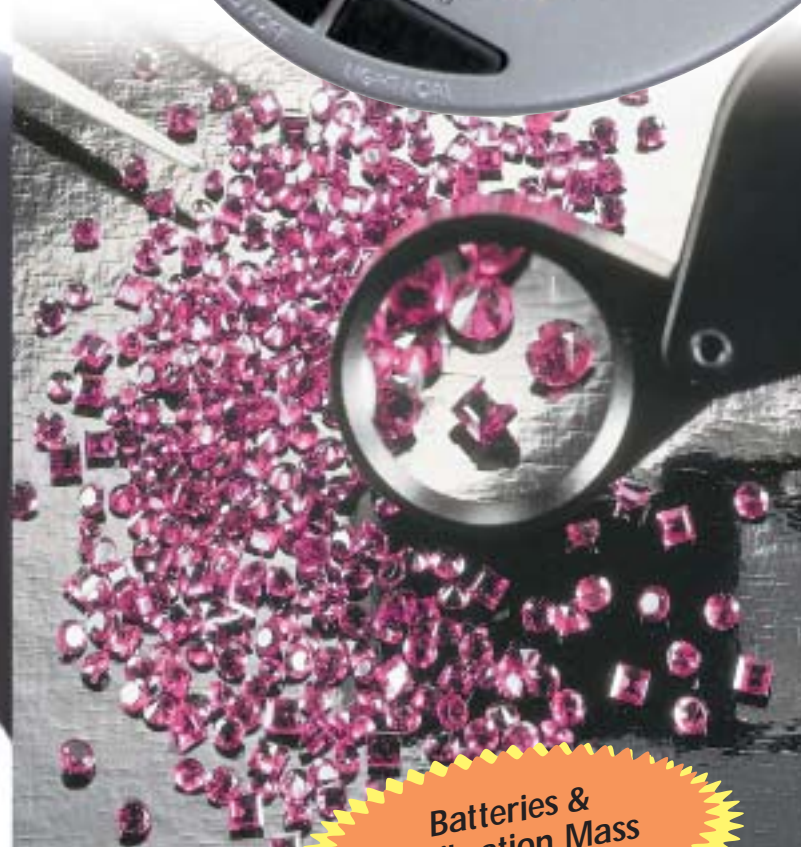

Batteries &  
Calibration Mass  
Included



**A&D** ...Clearly a Better Value  
A&D Company, Limited  
<http://www.aandd.jp>

# HJ-150

## ***/ Compact Design***

Easy to carry in your pocket, bag or briefcase and measurement can be done even in extremely cramped conditions (66 x 98mm)

## ***/ Light Weight***

Only 100g in weight (including batteries)

## ***/ Selectable Weighing Unit***

Meet any weighing application by choice of weighing unit of grams or ounces; simply change unit by pressing a button on the under-side of the scale

## ***/ High Weighing Resolution***

Precision weighing of 150 x 0.1g or 5.29 x 0.01 oz

## ***/ Tare Function***

One touch tare function enables the user to use a container for net weighing

## ***/ Battery Operation***

Uses 2 AAA alkaline batteries for easy transportation

## ***/ Auto Power Off Function***

To save battery consumption for a long battery life, it automatically turns itself off after 30 seconds of the display showing "0.0g", and after 2 minutes when the weighing is stable

## ***/ Large LCD Display***

Easy to read the 8mm high 7-segment LCD to prevent reading errors

## ***/ Backlight Function***

Beautiful blue backlight makes reading the display easy even in dark weighing conditions

## ***/ Standard Calibration Mass***

Standard calibration mass of 100g comes provided with the HJ-150 in its carton box to save you from purchasing calibration mass separately

## ***/ Calibration Function***

One-touch calibration with the standard calibration mass provided makes scale calibration extremely easy for accurate weighing

## ***/ Protective Cover***

Protects the weighing pan from overload and shock-load when it covers the scale; the protective cover can also double as a weighing container to be used to measure small objects

### *HJ-150 Specifications*

	HJ-150
Capacity x Resolution	150 x 0.1g / 5.29 x 0.01 oz
Display	7 segment LCD with Backlight (8mm height)
Power	2 x "AAA" size batteries
Pan size	60 (W) x 64 (D) mm / 2.36 (W) x 2.52 (D) inch.
Dimensions	66 (W) x 98 (D) x 20 (H) mm / 2.60 (W) x 3.86 (D) x 0.80 (H) in.
Weight	Approximately 100g / 0.22 lb



# AND

*...Clearly a Better Value*

### **A&D Company, Limited**

3-23-14 Higashi-Ikebukuro, Toshima-ku, Tokyo 170-0013 JAPAN  
Telephone:[81](3) 5391-6132 Fax:[81](3) 5391-6148  
<http://www.aandd.jp>

### **A&D ENGINEERING, INC.**

1555 McCandless Drive, Milpitas, CA. 95035 U.S.A.  
Telephone:[1](408) 263-5333 Fax:[1](408) 263-0119

### **A&D MERCURY PTY. LTD.**

32 Dew Street, Thebarton, South Australia 5031 AUSTRALIA  
Telephone:[61](8) 8301-8100 Fax:[61](8) 8352-7409

### **A&D INSTRUMENTS LTD.**

Unit 24/26 Blacklands Way Abingdon Business Park,  
Abingdon, Oxon OX14 1DY United Kingdom  
Telephone:[44](1235) 550420 Fax:[44](1235) 550485

### **<German Sales Office>**

Große Straße 13 b 22926 Ahrensburg GERMANY  
Telephone:[49](0) 4102 459230 Fax:[49](0) 4102 459231

### **A&D KOREA Limited**

Manhattan Bldg. 8F, 36-2 Yoido-dong, Youngdeungpo-gu, Seoul, KOREA  
Telephone:[82](2) 780-4101 Fax:[82](2) 782-4280